
XDCOLLECTION

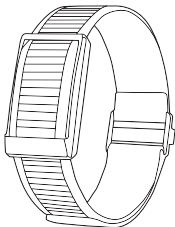Endo
Wearable health tracker

Copyright© XD P330.90XX

Contents

English	p. 3
Nederlands	p. 7
Deutsch	p. 11
Français	p. 15
Español	p. 19
Svenska	p. 23
Italiano	p. 27
Polski	p. 31
Explanation WEEE logo:	p. 35


ENGLISH

DEVICE REQUIREMENTS:

IOS 9.0 or later, Android 5.0 or later, Bluetooth 5.3.

INSTALL THE APP

Download the app by either scanning the QR code or search for "FitCloudPro" and download in the Apple Store or Google Play Store.

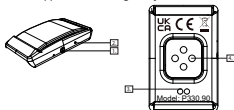DEVICE LAYOUT

1. On/Off
2. Indicator light
3. Charging port
4. Heart rate & blood oxygen sensor

SPECIFICATIONS

Battery Capacity	145 mAh
Input	5V/0.5A
Working time	12 days
Charging time	2 hour
BT version	5.4

BT pairing name	Endo
Compatibility	IOS 9.0 + Android 5.0
Waterproof	1ATM

*Times are approximate times only, depending on battery condition.

Important: Only use included accessories or certified and tested accessories to power this device. Using third party accessories can influence the item safety.

CHARGING

Before using the health tracker for the first time, please ensure it is fully charged. Use the equipped magnetic charger cable. Plug the USB-A end into a USB-A port on a computer or a compatible 5V/0.5A power adapter. While charging the light indicator is red. When fully charged the indicator light is off. Full charging takes approximately 2 hours.

CONNECTING

1. After you downloaded the app, turn on Bluetooth on your phone. Open the FitCloudPro app, enter the device page and press 'Bind Device'. The app will find the health tracker as 'Endo'.
2. After the health tracker is successfully paired with the app, a prompt will pop up asking for permissions to use the wireless call function. Press 'Allow'. Only then can the phone call function be used.

Note:

- Please do not directly pair the device in the Bluetooth list of the mobile phone system to avoid the situation where the "FitCloudPro" APP fails to connect to the health tracker.
- Please allow the "FitCloudPro" APP to obtain permissions such as "Location, Phone, Camera" of your mobile phone to ensure a successful pairing.

-
- Device must connect through the App. Do not pair directly by BT on your mobile device.

FUNCTION INTRODUCTION

1. Step Tracking

Tracks daily steps to monitor activity levels.

2. Heart Rate Monitoring

Automatically measures heart rate (approximately every 30 minutes).

3. Blood Oxygen Monitoring

Measures blood oxygen levels on demand.

4. Notifications

Vibrates and displays alerts for calls, messages, and apps.

5. Alarm Function

Set alarms for reminders and wake-up alerts.

6. Women's Health Tracking

Tracks menstrual cycles and related data.

7. Remote Camera Control

Shake the device to take photos on your phone.

8. Find Device

Locate your tracker via the app (device will vibrate).

9. Sport Modes

Supports activity modes (excluding water sports).

10. Sleep Tracking

Tracks sleep when worn overnight (recommended: 10:00 PM – 8:00 AM).

WARNING!

1. The app is not a medical device. Data obtained during the use of the app is for reference only, and should not be used for clinical diagnostic, medical research, diagnostic, or treatment purposes.
2. Please avoid the bracelet from strong impact, extreme heat and exposure.
3. Do not disassemble, repair or modify the watch yourself.
4. This watch uses 5V/0.5A for charging. It is strictly prohibited to use over-voltage load power supply for charging
5. Do not charge for more than 8 hours.
6. The operating environment of this watch is 0 ° - 40 °. It is forbidden to throw it into fire to avoid causing explosion.
7. Do not touch chemicals such as gasoline, cleaning solvents, propanol, alcohol, or insect repellent.
8. Do not use this product in high voltage and high magnetic environment.
9. If you have sensitive skin or you wear your bracelet too tightly, you may feel uncomfortable.
10. Please wipe off the sweat drops on your wrist in time. The strap may be exposed to soap, sweat, allergens or pollution components for a long time, which may cause itchy skin.
11. If you use it often, it is recommended to clean the bracelet every week. Wipe with a damp cloth and remove oil or dust with a mild soap. It is not advisable to wear a bracelet to take a hot bath.

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EXPLANATION WEEE LOGO:

ENGLISH

The crossed-out wheeled bin symbol indicates that the item should be disposed of separately from household waste. The item should be handed in for recycling in accordance with local environmental regulations for waste disposal. By separating a marked item from household waste, you will help reduce the volume of waste sent to incinerators or land-fill and minimize any potential negative impact on human health and the environment.

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EU Declaration of Conformity

ENGLISH

XD Connects B.V. hereby declares that this product is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU.

The complete Declaration of Conformity can be found at www.xdconnects.com > search on item number.

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